

R-Endo[®] : Protocole / Protocol

Vitesse 300 à 400 tr/min - Speed 300-400 rpm


1 Solvant
Solvent



2 Rm




3 Re

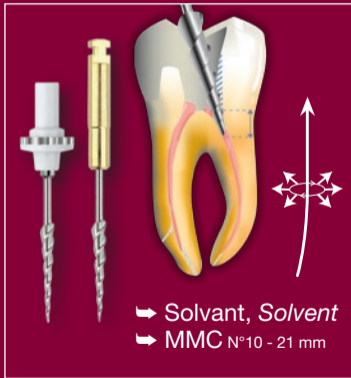

R-Endo[®] : Protocole / Protocol

Vitesse 300 à 400 tr/min - Speed 300-400 rpm

4 R1

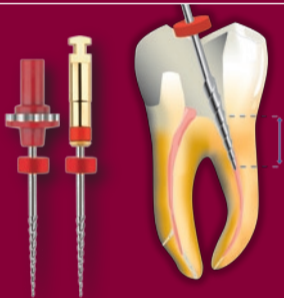


↪ NaOCl
↪ MMC N°10 - 21 mm
LT ? / WL ?

5 R2



↪ NaOCl
↪ MMC N°10 - 21 mm
à LT / to WL

6 R3

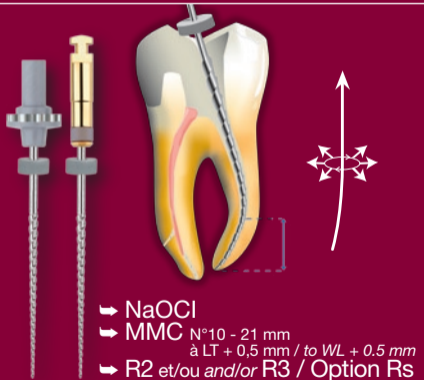

↪ NaOCl
↪ MMC N°10 - 21 mm
à LT + 0,5 mm / to WL + 0.5 mm
↪ R2 et/ou and/or R3 / Option Rs

